Q1: Missing parts when received.

A: If you find some parts are missing from the box, please kindly send us messages to let us know the name of missing part regarding to the first page of instructions. We will send replacement to you asap.

Q2: Some parts damaged when arrived.

A: Please send clear photos about damaged part to us as attachment, we will check and offer you best solution soon.

${\bf Q3:}\;\;$ The chair is stuck on lifting or reclining position, and no response when you press the lifting or reclining button.


A: Please kindly check if the power source(F) are connected well with lift motor. Sometimes the power source (F) may fall off from the motor.

If the connection is good, then check the lift motor if it has some anomaly, like some noise comes out when you press the lift button. If you can make sure the motor is defective, please contact us for replacement.

How to remove the chair?


Please note that you need to take apart the backrest before.

1. Force apart the buckles while take the backrest out from the arms.

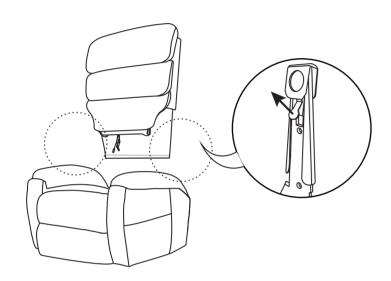


Please make sure you've got everyting lists as following before assembly.



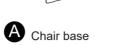




B Backrest

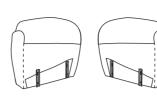

C Left and right armrest

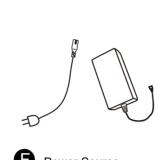

Power Source



Power cord


G Lift control



Remote control


Supporting Feet x4

Power Lift Chair Assembly Instruction

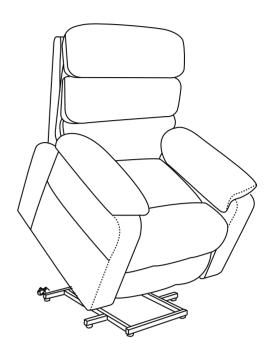


STEP 1

Find a small bag in the accessories box. Take out Supporting Feets from the small bag. Assemble the four Supporting Feet as shown on the bottom of mechanism $\!\scriptscriptstyle \circ$

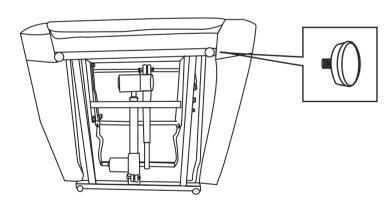


STEP 2


Insert two armrests(C) to the chair base(A), make sure the slots on the armrests fully assembled to the bottom part.

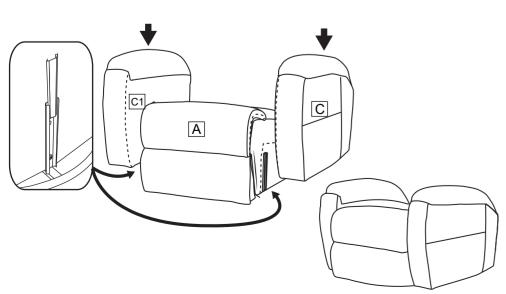



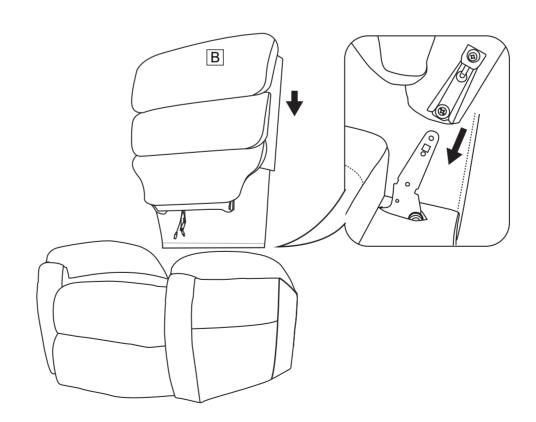

Insert chair back to the chair base, make sure the on the back fully assembled to the bottom part.

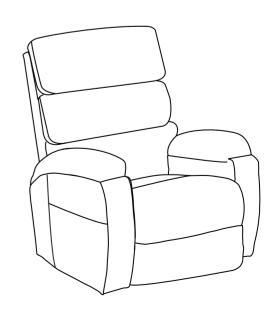


STEP 4

NOTE

- 1: Two same color wires connect.
- 2: Please connect all wires BEFORE plug power cord and power source into outlet.

No.1 and No.2 wire connect with the backrest and chair base.

- No.3 wire connects with backrest and power cord(D).
- No.4 wire connects with backrest and remote control(E). No.5 wire connect the motor and the wire from lift control(G).
- No.6 wire connect the motor and power source (F).

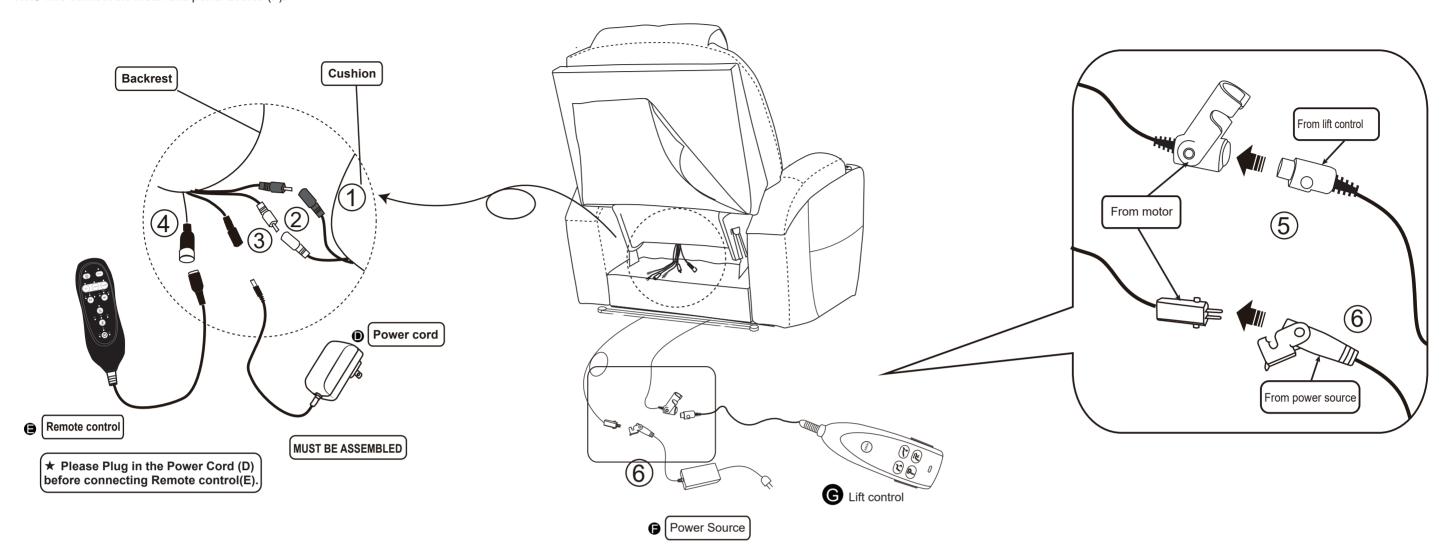

How to recline?

seat position

How to get footrest latch down and stand?

Turn on/off the massage function. You can select different massage parts to enjo You can select different massage mode to enjo Control the massage intensity from low to high. Press the button to rise the seat and Turn on/off the heating function. Press the button to descend the seat and w Z retrieve footrest. Select the heating temperature from M1 to M3. unfold the footrest. 15 **(1)** 40 -Press the button to Set a timer for massage from 15mins to 60mins. Press the button to lift the backrest. recline the backrest. 5 HOME button to return chair from any position to stand position.

stand position

How to use remote control?