

### RATED LOAD (lbs.)


1,000 lbs. DO NOT EXCEED RATED LOAD

### INSTALLATION

1. Read and follow all labels, instructions, and warnings prior to use. For additional safety information, or to obtain a copy of labels, instructions, and warnings, call 1-800-258-8503, or visit [www.ezaccess.com](http://www.ezaccess.com).
2. Use hand hold to transport and place Curb Plate; use caution to avoid injuring fingers when placing or removing Curb Plate.
  - ▲ Do not let feet or loads fall into, or otherwise become unstabilized by the hand hold.
3. Position top lip transition securely on support surface, ensure angle stops fully engage curb before and during use (FIG. 1).
4. Place the Curb Plate's top lip transition on upper support surface. Ensure that top lip transition fully overlaps the supporting surface and that the Curb Plate is perpendicular (90°) to the edge of the support surface. If obstructions prevent the top lip transition from fully overlapping or being placed at 90° from support surface, reposition Curb Plate as required.
5. Keep Curb Plate surface clean and free of debris. Clean oil or grease off of Curb Plate immediately with mild detergent and warm water (avoid alkaline cleaners).

### STORAGE

1. Curb Plate may be stored on its side or in any way that does not cause damage to the Curb Plate or property, or create a potentially dangerous situation.

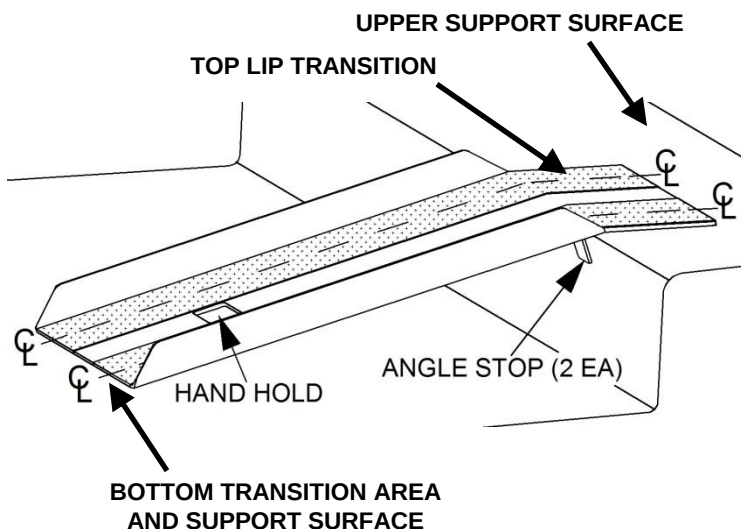

**FIG. 1**

### WARNING

- ▲ DO NOT EXCEED RATED LOAD.
- ▲ Do not drop Curb Plate.
- ▲ Do not use Curb Plate if angle stops are not fully engaged with support surface (FIG. 1).
- ▲ Avoid dragging Curb Plate as this will rapidly wear the Curb Plate and shorten useful life.
- ▲ Always exercise caution when using Curb Plate.
- ▲ Avoid working alone – always use plate with qualified assistants.
- ▲ Use only if both top lip transition and bottom transition area of Curb Plate are supported by level, stable surfaces.
- ▲ Use only if support surfaces are capable of supporting, at minimum, the same loads as the Curb Plate.
- ▲ Use only if sufficient maneuvering room is available at both ends of the Curb Plate to assure that all loads will enter, exit, and travel along the centerlines (Ⓢ) of the plate's skid-resistant tread. (FIG. 1).
- ▲ Do not let feet or loads fall into, or otherwise become unstabilized by the hand hold.
- ▲ Do not use as a vehicle ramp.
- ▲ Curb Plate will conduct heat rapidly. Do not place on or near open flame or hot objects as burns may result when using or handling.
- ▲ If Curb Plate is exposed to heat, hot weather or direct sunlight, use gloves when handling.
- ▲ Curb Plates in subfreezing temperatures will readily freeze any water that falls on the walking surface, thereby increasing the risk of slipping.
- ▲ Curb Plates reflect light. Use sunglasses or reposition Curb Plate to minimize unsafe glare during use.
- ▲ Use only if ice, snow, wet leaves and/or other debris, which may decrease traction and increase the risk of slipping, have been removed from Curb Plate.
- ▲ Water on the Curb Plate surface will decrease traction and increase risk of slipping. Use extreme caution when wet.
- ▲ When using Curb Plate in wet conditions, reduce weight, limit travel speed, and reduce Curb Plate angle if possible.
- ▲ If angle, load, or traction conditions are deemed unsafe, DO NOT USE Curb Plate!
- ▲ Inspect Curb Plate regularly for damage, including bending and cracked or broken welds as well as missing or ripped applied non-skid. If any part of the Curb Plate is unstable, damaged, or defective, DO NOT USE and contact your dealer.
- ▲ Aluminum Curb Plates will conduct electricity. Do not use or handle around live wiring or during electrical storms.
- ▲ Keep Curb Plate surface clean and free of debris.
- ▲ Clean Curb Plate with mild detergent and warm water (consult manufacturer before using abrasives or chemical cleaners).
- ▲ Clean any oil or grease immediately.

3-Year Limited Warranty. Please register at [www.ezaccess.com/warranty-satisfaction](http://www.ezaccess.com/warranty-satisfaction).

© EZ-ACCESS®, a division of Homecare Products, Inc. All rights reserved.  
All text and images contained in this document are proprietary and may not be shared, modified, distributed, reproduced, or reused without the express written permission of EZ-ACCESS.

